

Sartopore 2 0.45 µm

Bioburden & Particle Reductive Filter Cartridges & Mini Cartridges



Description

Sartopore 2 0.45 µm rated filter cartridges are designed for bioburden reduction and particle removal from a broad range of pharmaceutical products. They offer extremely high flow rates and total throughputs and are therefore ideally suited for membrane prefiltration of aqueous solutions and highly viscous, difficult to filter pharmaceutical products.

Applications

Typical applications include bioburden reduction and particle removal from:

- Buffers
- Biological Fluids
- Ophthalmics
- LVP
- Antibiotics
- Bulk pharmaceutical products

Compatibility

Featuring a unique hydrophilic polyethersulfone membrane, Sartopore 2 0.45 µm cartridges are compatible with solutions from pH 1 to pH 14 and are unaffected by numerous steam sterilization cycles. They are ideally suited for filtration of solutions with high | low pH and for multiple SIP | CIP cycles.

Performance

Sartopore 2 0.45 µm cartridges provide exceptional high flow rates, resulting in economical sizing of filtration systems. Due to the "built-in prefiltration" by a 0.8 µm membrane, Sartopore 2 0.45 µm rated cartridges offer outstanding total throughputs.

Wettability

Sartopore 2 cartridges can be easily wetted out for integrity testing even after drying at 80 °C for 12 hours.

Microbiological Retention

Sartopore 2 0.45 µm rated filter cartridges are validated for removal of *Serratia marcescens* with a Log Reduction Value (LRV) of 7 according to HIMA and ASTM F-838-83 guidelines.

Quality Control

Each individual element is integrity tested by diffusion and bubble point test prior to release, assuring absolute reliability.

Documentation

Sartopore 2 cartridges are designed, developed and manufactured in accordance with an ISO 9001 certified Quality Management System. A Validation Guide and Extractables Guide are available for compliance with regulatory requirements.

Specifications

Materials

Prefilter Membrane:	Polyethersulfone, asymmetric
Endfilter Membrane:	Polyethersulfone, asymmetric
Support Fleece:	Polypropylene
Core:	Polypropylene
End Caps:	Polypropylene
O-Rings:	Silicone (optional EPDM or Viton)

Pore Size

0.8 µm + 0.45 µm

Available Sizes | Filtration Area

Standard Cartridges

Size 1	10"	0.6 m ² 6.5 ft ²
Size 2	20"	1.2 m ² 12.9 ft ²
Size 3	30"	1.8 m ² 19.4 ft ²

Mini Cartridges

Size 7	0.05 m ² .5 ft ²
Size 8	0.1 m ² 1 ft ²
Size 9	0.2 m ² 2 ft ²

Available Adapters

21, 25, 27, 28

Operating Parameters

Max. Allowable Differential Pressure:	5 bar 75 psi at 20 °C
	2 bar 29 psi at 80 °C
Max. Allowable Back Pressure:	2 bar 29 psi at 20 °C

Integrity Test Limits

Maximum Allowable Diffusion at 1.7 bar / 22 psi & 20 °C

Cartridge Size	Maximum Diffusion	Minimum Bubble Point
Size 1	12 ml/min	2.2 bar/32 psi
Size 2	24 ml/min	2.2 bar/32 psi
Size 3	36 ml/min	2.2 bar/32 psi
Size 7	3 ml/min	2.2 bar/32 psi
Size 8	4 ml/min	2.2 bar/32 psi
Size 9	6 ml/min	2.2 bar/32 psi

Membrane Filter Cartridges

Specifications

Extractables

Sartopore 2 0.45 µm rated filter cartridges meet, or exceed the requirements for WFI quality standards set by the current USP.

Regulatory Compliance

100% Individually integrity tested

Integrity test correlated to HIMA/ASTM F 838-83 Bacteria Challenge Test using *Serratia marcescens*

Non-pyrogenic according to USP Bacterial Endotoxins

Meets USP Plastics Class VI biological reactivity test, in vivo

Non-fiber releasing according to 21 CFR

Sterilization

In-Line Steam Sterilization:

134 °C, 20 min. at max differential pressure of 0.5 bar | 7.25 psi

Autoclaving:

134 °C, 2 bar | 29 psi, 30 min

Sterilization Cycles

In-Line Sterilization: Min. 25
Autoclaving: Min. 25

Technical References


Validation Guide:

SPK 5732-e

Extractables Guide:

SPK 5731-e

Ordering Guide:

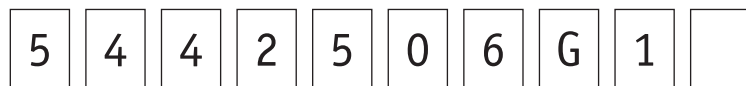

Filter Family
544 = Sartopore 2


Filter Format
21 = Double Open End
25 = Standard Cartridge (Code 7)
27 = Standard Cartridge (Code 3)
15 = Mini Cartridge

Final Membrane Pore Size
06 = 0.45 µm

Pre-filter Membrane Pore Size
G = 0.8 µm


Cartridge Size | Package Quantity
1 = 10" Standard Cartridge, 0.6m²/6ft², Single Unit
2 = 20" Standard Cartridge, 1.2m²/12ft², Single Unit
3 = 30" Standard Cartridge, 1.8m²/18ft², Single Unit
7B = 500 cm² / .5ft² Mini Cartridge, Pack of 5
8B = 0.1 m² / 1ft² Mini Cartridge, Pack of 5
9B = 0.2 m² / 2ft² Mini Cartridge, Pack of 5

Water Flow Rates for Standard Cartridges

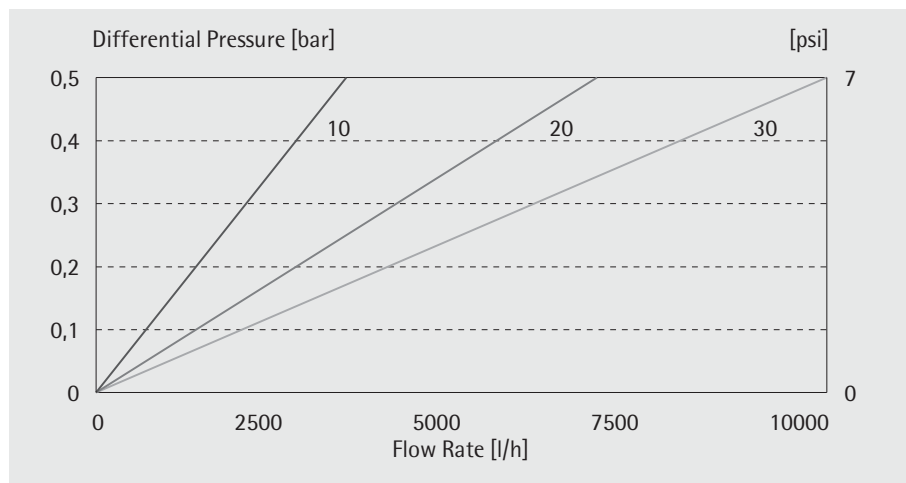

Standardized at 20 °C

Total Throughput Comparison

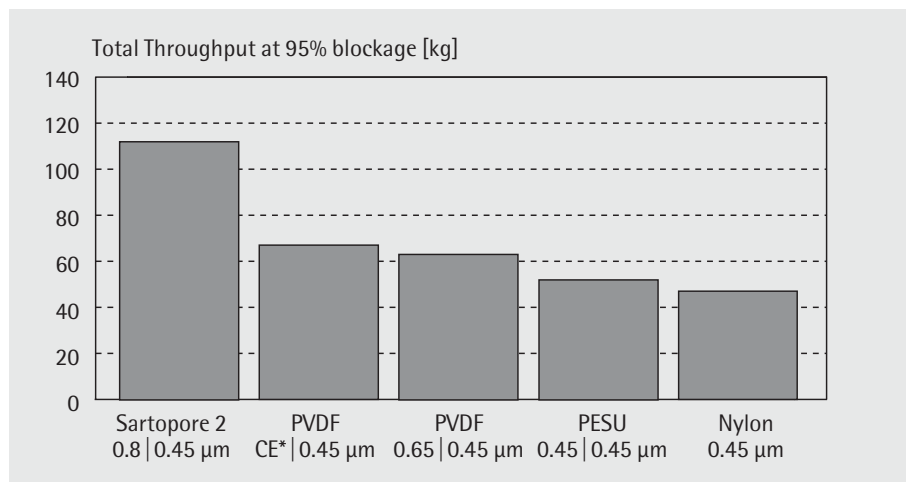

10" Cartridges

* Cellulose Ester prefilter

Sartorius AG
Weender Landstrasse 94-108
37075 Goettingen, Germany
Phone +49.551.308.0
Fax +49.551.308.3289
www.sartorius.com

USA +1.631.2544249
UK +44.1372.737100
France +33.1.69192100
Italy +39.055.634041
Spain +34.91.3586100
Japan +81.3.33293366

Specifications subject to change without notice.
Printed and copyrighted by Sartorius - W/sart-000 - G
Publication No.: SPK2039-e04041
Order No.: 85030-513-44